



# Medical cart compact-cart

NEW COLORS



## TECHNICAL DETAILS

Load capacity	Basic frame "Economy"	Basic frame "Profi"
Total load capacity	150 kg / 330 lbs	180 kg / 396 lbs
Load capacity of shelf	50 kg / 110 lbs distributed load	50 kg / 110 lbs distributed load
<b>Dimensions - Basic frame</b>		
Width in mm	555 / 605	618 / 668
Depth in mm	625	624
Total height in mm	842 / 1130 / 1290 / 1450	1146 / 1306 / 1466
<b>Dimensions - Shelf</b>		
Width in mm	450 / 500	450 / 500
Depth in mm	491	491
<b>Electrification</b>		
Isolating transformer	no	1200 VA / 1600 VA / 2000 VA
Central on / off switch	no	yes
Socket strips	yes	yes

Get in touch >>

Europe: +49 89 614425-0 | sales@itd-cart.com  
 America: +1 800 947-3901 | salesusa@itd-cart.com  
 China: +86 21 5031 8818 | saleschina@itd-cart.com  
 India: +91 86 0015 4064 | salesindia@itd-cart.com  
 Australia: +61 3 9330 5050 | salesaustralia@itd-cart.com

www.itd-cart.com



Basic frame "Profi"




DS-0288.799 - 10/2024. Printed in Germany. All rights to make changes reserved.

## CUSTOMIZED DESIGN

### New colors

#### BASE COLOR



RAL 9016  
Traffic white

#### DECORATIVE TRIM COLORS



RAL Design 2507025  
Meadow blue



RAL 5017  
Traffic blue



RAL 7016  
Anthracite grey

### Color design

For an appealing color design, decorative elements such as handles and corners can be powder-coated in selected colors. It is also possible to implement your corporate design color.



# HIGH PAYLOAD & DIFFERENT BASE FRAMES

**Different base frames –**  
in different heights and widths with a payload of up to 180 kg / 396 lbs




Basic frame "Economy"



Basic frame "Profi"

# HIGH STABILITY FOR LARGE MONITORS

**Monitor holders –**  
for monitors up to 18 kg / 39.6 lbs

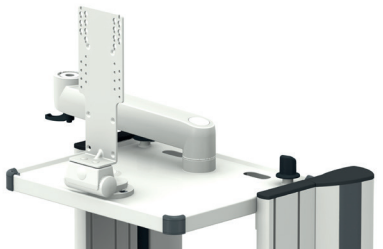

Monitor shelf, double for flexion-port support arm



VESA 75 / 100, tilt and swivel





VESA 75 / 100 adjustable, tilt and swivel




**Air-permeable rear door –**  
prevents heat congestion of the equipment and unwanted access

# SAFE ELECTRIFICATION

**Separate data and power cables –**  
prevent electromagnetic interference

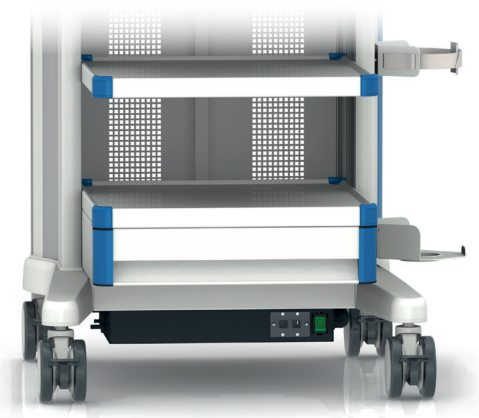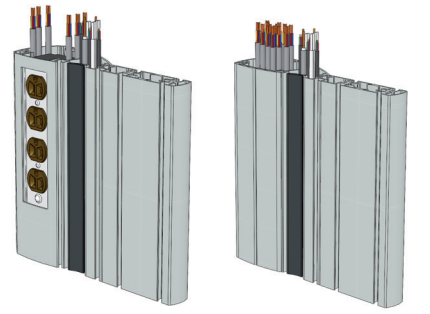

**Isolating transformer up to 2000 VA –**  
electrical safety for patients and users, international use by switching the voltage

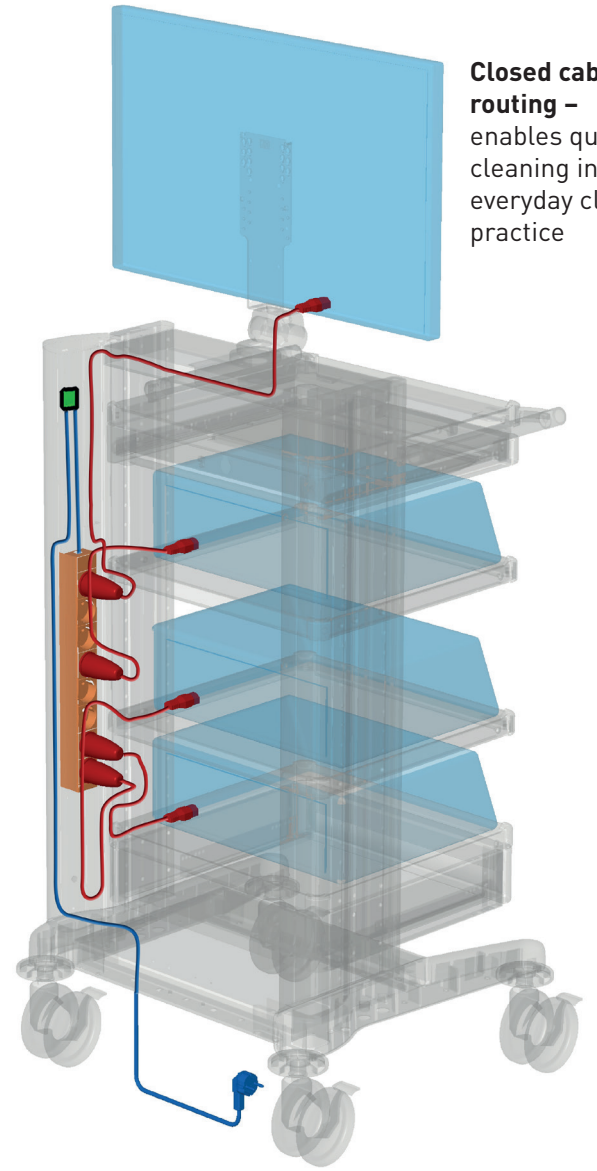

**Closed cable routing –**  
enables quick cleaning in everyday clinical practice